


# Smile CHANGERS

ADVANCED TEETH WHITENING SYSTEM

## USER MANUAL

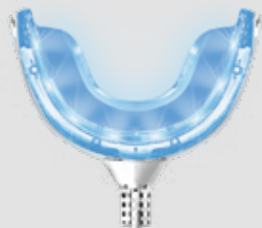

WWW.SMILECHANGERS.NET

IPHONE | ANDROID COMPATIBLE

Whiten your teeth in the comfort of your own home with Smile Changers Advanced Teeth Whitening System. Smile Changers smart teeth whitening mouth tray LED light adopts special cold light technology, combined with mouthpiece and high intensity cold light. Blue ray reacts with the professional teeth whitening gel, which can remove the stains adhered on the teeth's surface.

The soft mouthpiece is designed based on the principles of human oral engineering. You will feel more comfortable in the course of teeth whitening process. When you bite the mouthpiece, it will completely fit and cover fit and cover your teeth.

Recommended by #1 Cosmetic Dentist, Dr. Heavenly.

### CERTIFICATION



FDA



CE




CPSR




RoHS



Rinse mouth and clean teeth

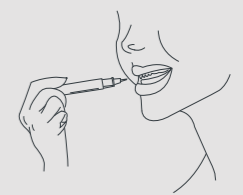

Apply whitening gel to teeth's surface

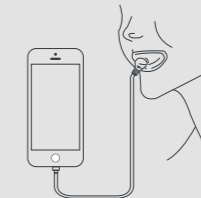

Connect Apple/Andriod phone or USB port



16 minutes treatment



Rinse your mouth



Clean mouth tray

# Smile CHANGERS

ADVANCED TEETH WHITENING SYSTEM

### STEPS FOR USAGE

1. Please rinse tray before using, then use teeth wipe to dry the teeth.
2. Look in mirror on the back of power bank, comparing your teeth with shade guide paper, record which stage you are in.
3. Twist teeth whitening pen until the gel flows to the brush. Apply teeth whitening gel evenly on your teeth (thickness approximately 1mm). Avoid applying teeth whitening gel directly on your gum or inside of your lips!
4. Connect the smart teeth whitening device with power bank, put mouthpiece into your mouth and bite tightly when the light is on. Take deep breaths before treatment to decrease saliva production. (Excessive amount of saliva will dilute the gel's concentration.)
5. Take out mouthpiece after 16 minutes. Rinse your teeth with wan water and compare your teeth with shade guide paper again to see how many stages you have improved.
6. The smart teeth whitening device can be put directly in warm water or washed under running water.

**INGREDIENTS:** Medicinal Hydrogen Peroxide, Carbomer, SiO<sub>2</sub>, WDJ-38 Glycerinum, Potassium Sorbate, Xylitol, Deionized Water

### SUITABLE RANGES

Yellow teeth or badly-stained teeth caused by smoking, drinking tea or coffee, and other reasons.

Age-related yellow teeth or badly-stained teeth.

Light or moderate tetracycline pigmentation teeth and dental fluorosis.

Those who want a whiter and brighter smile.

**STORAGE:** Avoid sunlight and keep in cool place.

**SHELF LIFE:** 12 months

### WARM PROMPT

1. Rinsing thoroughly, water on the teeth's surface should be dried entirely, applying gel to teeth evenly.
2. Excessive saliva will dilute the gel's concentration.
3. Everyone's teeth are different. Few people will feel slight sensitivity, which is normal phenomenon. You can apply whitening repair essence to ease discomfort. If your sensitivity is severe, please reduce the use time or stop using immediately.

4. Avoid directly applying teeth whitening gel to gum area - wipe off excess.
5. If gel is on your face or clothes, you should wash them without any delay.
6. Please put cap on the top of syringe after finishing teeth whitening treatment, which can avoid active ingredient's volatilization.
7. Please store our product in a place which is out of children's reach.

### CAUTION

Not suitable for ceramic teeth and false teeth.

Not suitable for teeth discoloration caused by lesion or medicine.

Not suitable for severe tetracycline and decayed teeth.

Not suitable for defective enamel, exerted dentin and damaged teeth.

Not suitable for children under 12 years old and pregnant women.